

Wash Bottles

Clear and informative product safety information

Wash bottles safety labeled self venting

The information provided on the front of each bottle is specific to the product contained and includes the following — the Chemical Abstract Service (CAS) No., National Fire Protection Association (NFPA) codes and symbols, Protective Equipment Symbol (PES), Target Organ Information and Hazard Codes.

You also have the option of using other liquids within the bottles and creating the appropriate safety information on the outside of the bottle by using the 'self labeling' versions.

The closure and stem are molded in one piece to make them leak-proof. The angled neck and tip provide a uniform stream of liquid and the tip can be cut back if required to increase the flow

- Wide bottle mouth makes filling easier
- Non-metallic vent that prevents the risk of solvent leakage
- With the exception of the bottle designed for Sodium Hypochlorite, which is a natural LDPE bottle for light protection and has a HDPE closure, the bottles are made from LDPE, with a PP closure
- Six named chemicals, a self labeling** and assorted pack of bottles to choose from
- 500 mL bottle size, supplied in packs of 6

| Part No. | Description | Color |
|-----------------|---------------------|--|
| HS120250 | Acetone | Red <input checked="" type="checkbox"/> |
| HS120251 | Distilled Water | Natural <input type="checkbox"/> |
| HS120252 | Isopropanol | Yellow <input type="checkbox"/> |
| HS120253 | Ethanol | Natural <input type="checkbox"/> |
| HS120254 | Methanol | Green <input type="checkbox"/> |
| HS120255 | Sodium Hypochlorite | Blue <input type="checkbox"/> |
| HS120256 | Self Labeling** | Natural <input type="checkbox"/> |
| HS120257 | Assorted* | <input checked="" type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> <input type="checkbox"/> * |

*Assorted pack includes one of each, acetone, distilled water, isopropanol, ethanol, methanol, sodium hypochlorite and a self labeling bottle

**Self labeling bottles must only be used with chemicals that are safe and compatible with LDPE





Wash bottles color-coded

Color-coded for easy identification

The closure and stem are molded in one piece to make them leak-proof.

- The angled neck and tip provide a uniform stream of liquid and the tip can be cut back if required to increase the flow
- Wide bottle mouth makes filling easier
- Color-coding for user, application or location identification
- Bottle is made from LDPE, with a PP closure
- 500 mL bottle size, supplied in packs of 6

| Part No. | Color | |
|-----------------|--------|--------------------------|
| HS120245 | White | <input type="checkbox"/> |
| HS120246 | Red | <input type="checkbox"/> |
| HS120247 | Yellow | <input type="checkbox"/> |
| HS120248 | Green | <input type="checkbox"/> |
| HS120249 | Blue | <input type="checkbox"/> |

Eco-smartFunnel™ Paper Funnel

A clean transfer every time


The Eco-smartFunnel™ is supplied flat-packed to minimize storage space and help keep areas clean. When ready to use, apply a minimal amount of finger pressure and the funnel takes its working shape. The Eco-smartFunnel™ is an environmentally friendly paper product that will naturally biodegrade in the standard waste stream.

- Fits a wide variety of openings and can be used with virtually any liquid or powder
- Innovative “backrest” allows hands-free pouring
- Maximum flexibility - cut the stem to any size or shape to meet with your needs
- Supplied in packs of 25
- Not recommended for hot liquids
- Not recommended for cleaning in autoclave
- Measures 194 x 140 x 51 mm



Innovative
‘backrest’ allows
hands-free pouring

| Part No. | Description | Pkg. |
|-----------------|--|-------|
| HS120245 | Eco-smartFunnel™ 186 mL cap. Top Dia. 6cm, Bottom Dia. 0.64 cm, Height | 25/pk |

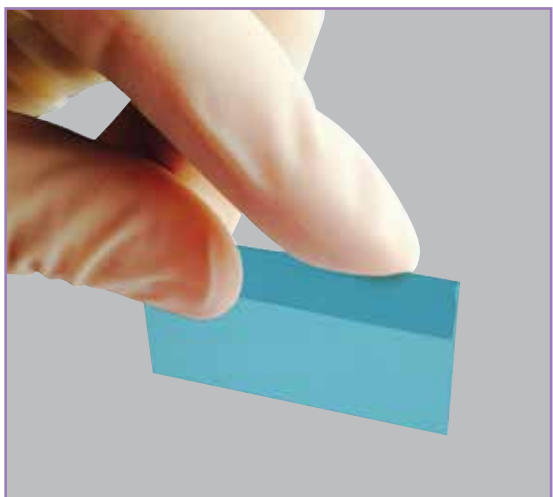

smart Slicer™ Laboratory Safety Razor Blade

Say No to dangerous sharps!

At just 0.22 mm the Patent Pending product is as thin as a single edge industrial razor blade, but much safer.

- Cleanly slice agarose gels right on the light box without scratching the surface
- Ideal for cutting soft tissue and plant matter
- Adjust the size required by simply cutting with scissors
- Available in two sizes and made in the USA from food grade PET
- Supplied in packs of 100

| Part No. | Size | Pkg. |
|-----------------|-----------------------------------|--------|
| HS120259 | 19 mm H x 38 mm W x 0.22 mm thick | 100/pk |
| HS120260 | 25 mm H x 50 mm W x 0.22 mm thick | 100/pk |