

VEE GEE®

Buret KIT

Three Great Buret Products all in one Kit.

- Two each of either: 10, 25 or 50mL Class A Burets
- One: Buret Retort Stand
- One: Double Buret Clamp with White Finish

VEE GEE 20114-Series Burets meet ASTM E-287 Class A, w/ PTFE Stopcocks.

Made from precision, uniform Schott DURAN® true-bore tubing to ensure the stipulated capacity tolerances. These burets feature durable, permanent markings, fine, sharp lines, and large, easy-to-read numbers. Supplied with PTFE stopcocks to eliminate freezing and lubricant contamination and to provide a precision fit within the barrel. All sizes are supplied with a 2mm stopcock bore plug.

Non-corrosive, molded Retort Stand provides an excellent replacement for metal ones. Even the metallic rod has a polypropylene covering making it non-corrosive.


The rod can be threaded into the heavy base giving way to a perfectly stable retort stand. Interfaces with the VEE GEE Class A Buret Series and Buret Clamp.

The VEE GEE white finish double buret clamp is made of aluminum, die casted and painted white.

Black Rubber sleeves allow for firm gripping of burets. The clamps are adjustable and hold two burets from micro to 100mL sizes.

The thumb style screw key attaches easily to support rods.

⚠ WARNING: This product can expose you to chemicals including lead, which is known to the state of California to cause cancer and birth defects or other reproductive harm. For more information go to: www.P65Warnings.ca.gov

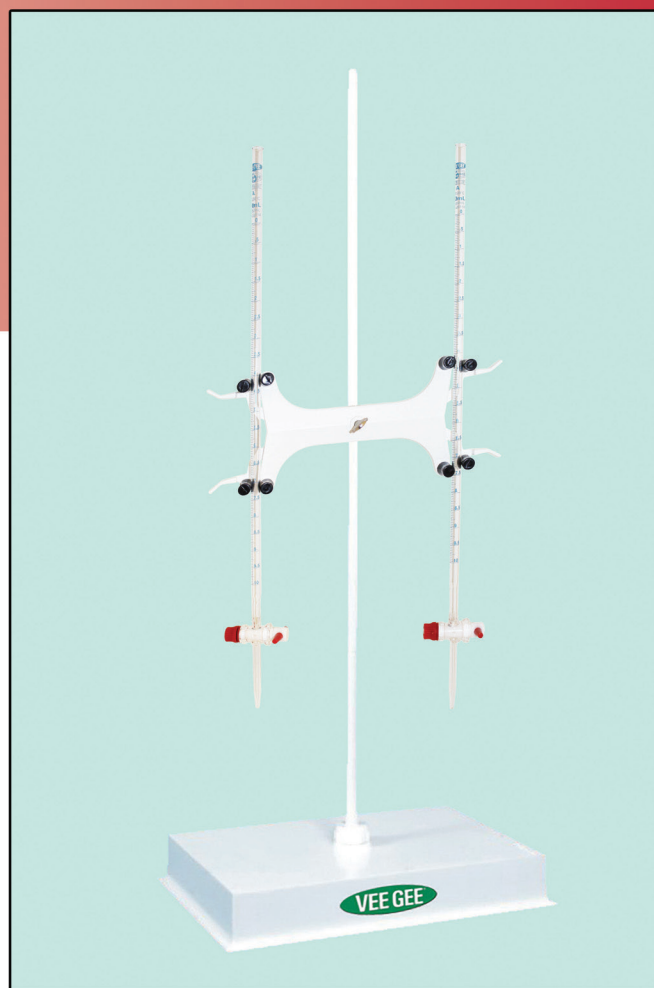

Cat. No.	Cat. No. Kit Components
2011BK-10	Two 20114-10 (10mL Burets), One 20117 (Retort Stand), & One 3021 (Clamp)
2011BK-25	Two 20114-25, (25mL Burets), One 20117 (Retort Stand), & One 3021 (Clamp)
2011BK-50	Two 20114-50, (50mL Burets), One 20117 (Retort Stand), & One 3021 (Clamp)



VEE GEE Scientific, Inc.